

Highlights

- Center Wheel Drive for maximum maneuverability
- Four wheel drive: freedom to explore new places
- Extreme safety and stability that only 4W MAX drive can deliver
- Best in class climbing up to 13 cm (5.1")
- Power adjustable front wheel height
- Powered Variable Ground Clearance "VGC"
- Multi axle suspension for maximum comfort on any surface
- Best in class Sedeo Pro modular seating system standard with total configurability and center point of gravity adjustment
- Proven quality components used throughout MAX
- High torque 4 pole motors
- Easy maintenance and service




To see all MAX features, more pictures and movies, visit our website or scan the QR code

Technical Data

Model: MAX 300	Metric	Imperial
Seat width (SW)	39-54	24.2 - 33.5
Seat depth	42-60	16.5 - 37.3
Seat height (min.)	43	26.7
Total width (W)	61	24
Total length (L)	104	41.3
Total height (incl. Backrest) (H)	106	41.7
User weight (maximum)	160	352
Speed	6, 10, 12	3.7, 6.2, 7.4
Battery capacity (Ah)	75 AGM, 73 Gel	75 AGM, 73 Gel
Range (with 73 Ah batteries)	32	19.9
Turning circle (D)	113	44.5
Dynamic stability / max. safe slope	15° (33%)	15° (33%)
Obstacle height	13	5.1

Standard configuration
Powered 4WD wheelchair
Choice of five standard colors
Powered Variable Ground Clearance
Powered front wheel height adjustment
Seat height 48 cm (18.9")
Speed 10 km/h (6.2 Mph)
Drive and front wheels air
Double rear castors puncture proof
Multi-axle suspension
PG Drives Technology R-Net controller
Sedeo Pro "quick adjust" seat frame
Profiled seat cushion
Padded backrest
Central leg rest
Folding backrest
Flip up armrests
Manually adjustable seat tilt 0,3 and 6°
Maximum user weight 160 kg
75 Ah AGM batteries
24V/8A charger
Vehicle tie-downs

Options / accessories
Headlights; Turn signals; Brake lights
Powered seat tilt up to 50°
Powered seat elevation up to 30 cm (11.8")
Powered backrest recline up to 160°
A large variety of backrests and seat cushions
Angle adjustable leg rests
High torque motors
2 colors + 190 RAL colors on request
A large range of special controls
....and many more (see the order form)

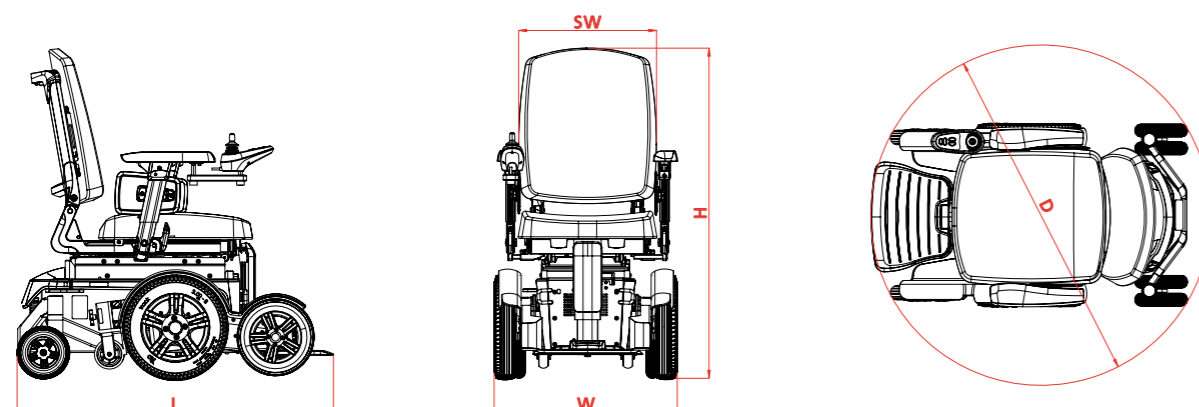

Distributor:

Manufactured by:

LEVO AG
Switzerland

Anglikerstrasse 20
CH-5610 Wohlen

Tel. +41 (0)56 618 44 11
Fax +41 (0)56 618 44 10

office@levo.ch
www.levo.ch

LEVO.dt
A Dane Technologies Company

LEVO, partners for a better life

LEVO.dt
A Dane Technologies Company

Discover the freedom, power, and comfort of MAX and four wheel drive

MAX 300: the "Go Anywhere" power chair that opens up a brand new world of opportunities for you to explore

MAX is the revolutionary new power chair from LEVO. It combines legendary mid wheel maneuverability with the power and safety of full time four wheel drive. MAX's innovative Swiss design is matched with the award winning Sedeo seating system that adjusts to your every need and desire. This powerful combination delivers the comfort and performance that easily conquers obstacles and opens up a whole new world of mobility opportunities.



Sedeo seat system with VR2 or R-net controls



Multi axle suspension system

Let MAX take you anywhere you want to go
With the MAX 300 and its powered Variable Ground Clearance "VGC", you have a partner that instantly adjusts to every situation. Its patent pending power front drive wheels raise or lower, allowing you to easily go up or down over curbs, up to 13 centimeters with maximum safety and stability. MAX's best in class performance is enhanced by its sophisticated suspension system that delivers maximum comfort and traction in every situation. All of this is combined with proven and dependable high torque 4 pole motors that power you up to 12 KM/H. The MAX delivers the maneuverability you need indoors with the safety, power, and comfort to explore your outdoor world.

combines a highly adjustable seat frame with a wide choice of support and positioning options, including a full range of high-quality seat cushions and backrests. Sedeo can easily be customized to meet your needs and requirements, delivering maximum comfort and support all day every day. With available options, tilt recline, elevation, elevating leg- and foot rest, you can always find the perfect position.

Best in class seat system
The MAX 300 is equipped with the best in class seat system: Sedeo Pro. This award winning seat system

Individualize MAX to your personal taste
When you choose MAX, we want you to have the opportunity to make your power chair reflect the unique individual you are. Choose any of the five vibrant colors that come standard or indulge, and choose one of our many custom colors. Just contact us and order your favourite RAL colour. With the MAX 300, the choice is yours.



Curbs are easy with MAX 300

